CARRERA

DESIGNED IN GERMANY

EN

WATER FLOSSER

Instruction manual

RU

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www.carrera.de

Welcome to the **CARRERA** club!

The story of **CARRERA** starts in the 80' in Germany as a lifestyle across two product categories: personal care and home appliances.

CARRERA engineers and designers were inspired by the legendary Carrera Panamericana, the longest and the most challenging race in the 50'. Only the racers who were ready to go beyond the limitations of their own bodies and materials stood a chance of winning.

Today, that same philosophy drives **CARRERA** to develop unique and powerful products.

Each product is made of hi-tech components and quality materials outlining power, quality, innovative features and effective design.

This approach resonated with millions of people, boosting **CARRERA** brand recognition in Europe up to 90 %. **CARRERA** products have become an embodiment of German precision and craftsmanship.

Today, **CARRERA** continues to invest into research and development. Each **CARRERA** product provides users with a sense of luxury, comfort and efficacy.

CARRERA is not just about tools. It is about a lifestyle.

CONTENT

01 SAFETY PRECAUTIONS	4
<hr/>	
02 INTENDED USE	5
<hr/>	
03 APPLIANCE OVERVIEW	5
<hr/>	
04 BEFORE USE	6
<hr/>	
05 OPERATION	8
<hr/>	
06 CLEANING & MAINTENANCE	9
<hr/>	
07 TROUBLESHOOTING	10
<hr/>	
08 DISPOSAL	11
<hr/>	
09 TRANSPORT AND STORAGE	11
<hr/>	
10 TECHNICAL SPECIFICATION	12
<hr/>	
11 SCOPE OF DELIVERY	12
<hr/>	
12 WARRANTY	12
<hr/>	
13 ADDITIONAL INFORMATION	13
<hr/>	

01 SAFETY PRECAUTIONS

Read this instruction manual carefully before using the appliance. Keep this instruction manual, the packing box and, if possible, other elements of the appliance's packaging. The manufacturer shall not be liable for any damage or malfunction of the appliance or for any damage resulting from incorrect or improper use of the appliance.


- Check the appliance and USB cable regularly for damage. Do not use the product if any damage is detected.
- Do not use the product for other than its intended use.
- Do not use the product outdoors.
- Do not immerse the product or USB cable in water or other liquids.
- Do not expose the product or USB cable to high temperatures, direct sunlight, dampness, or other adverse conditions.
- If you have recently undergone oral surgery or are undergoing oral treatment, it is not recommended to use the appliance without prior consultation with a doctor.
- If your gums start bleeding during use, you should stop using the product until you consult your doctor.
- Remove jewelry (rhinestones, grills, etc.) from the teeth before using the appliance.
- Always turn the appliance off before changing the nozzle or filling the water tank.
- Do not turn the appliance on when the water tank is empty. Do not fill the water tank with hot water.
- Do not leave the appliance running unattended, turn it off immediately after use.
- Be sure to pour any remaining water out after use. Clean the appliance after each use.
- Do not direct the water stream into the eyes, under the tongue, into the ear, nose or other sensitive areas.
- Do not use accessories that are not included in the scope of delivery.
- Do not subject the appliance to mechanical impact. Do not apply excessive pressure to the nozzle. Do not use deformed or mechanically damaged nozzles.
- Do not allow children to play with the appliance, the USB cable or the packaging.
- The appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience or knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety.

- Do not attempt to repair the appliance or replace any parts yourself. In case of malfunctions, contact the nearest service center.

02 INTENDED USE

The water flosser is designed for oral care and cleaning of the interdental space. The appliance is intended for use in commercial or domestic areas in accordance with the intended use and the safety precautions described in these operating instructions.

03 APPLIANCE OVERVIEW

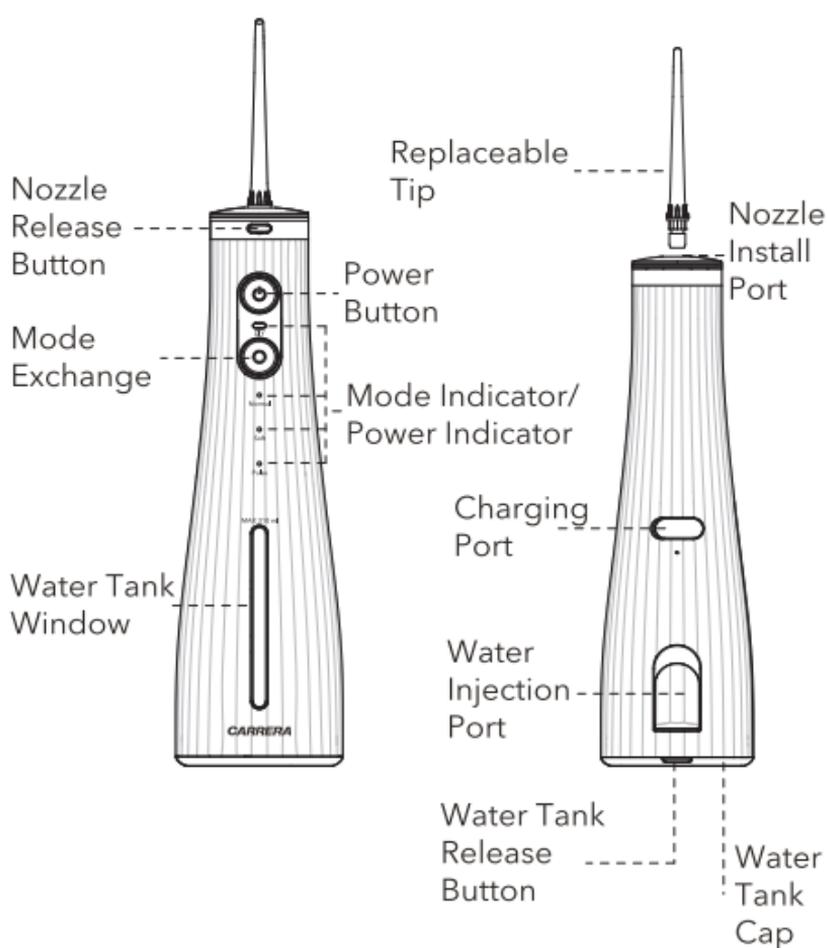04 BEFORE USE

Charging the appliance


- Charge the appliance for 4 hours before using it for the first time. This charging is necessary every time the appliance has not been used for a long time.
- Connect the appliance to a power supply (not included) via USB cable. The output power of the power supply should be no more than 15 W. The ambient temperature during charging should be 5-35 °C, otherwise the battery performance may deteriorate.
- You can determine the charge level by the mode indicators: when the battery level is low, all indicators will blink.
- During charging, the LEDs will blink slowly and will light up solid when the battery is fully charged.
- One full battery charge is enough for 35 days of use.
- After charging, disconnect the appliance from the power source immediately. The appliance cannot be used during charging.
- When the battery is fully charged, all LED-lights on will be steady on when the appliance is turned off. As the battery is used, the lights will gradually go out.

WARNING!


- **Do not leave the appliance discharged for long periods of time.**
- **When charging for the first time or after a long period of non-use, the indicators may not light up immediately after charging has started. Wait about 5 minutes, this is not a malfunction.**
- **Do not touch the USB cable with wet hands.**
- **If the appliance or USB cable become very hot during charging, disconnect the appliance from the power source immediately and contact an authorized service center.**
- **The appliance may vibrate slightly during charging, this is not a malfunction.**

Choosing and installing the tip

- The set includes 3 replaceable tips:

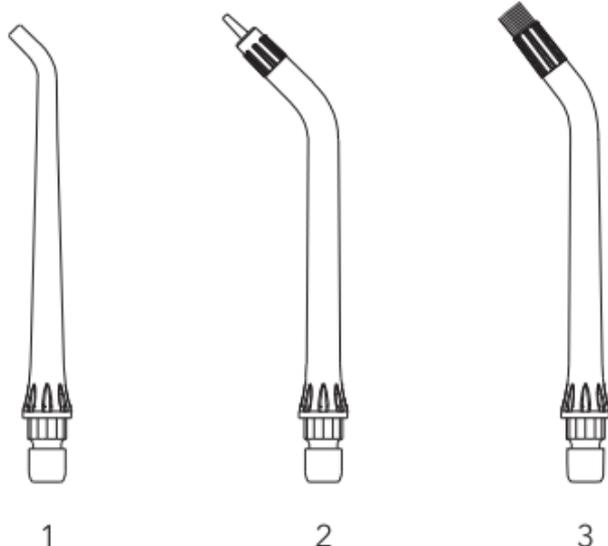
1. Standard tip for daily oral care.
 2. Periodontal tip for inflamed gums (suitable for periodontitis, periodontal disease and other gum diseases).
 3. Orthodontic tip for cleaning braces as well as bridges, implants and hard-to-reach areas of the mouth.
- Insert the tip into the special hole in the top of the housing until it clicks into place. Make sure it is securely locked in place.
 - Press the button on the body of the appliance to remove the tip.



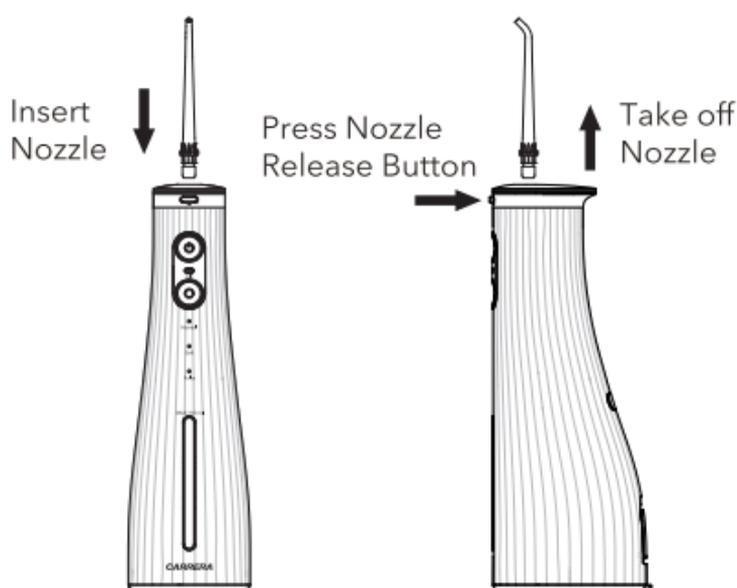
Filling the water tank

Fill the tank with clean water through the opening in the body. Close the cover of the opening securely.

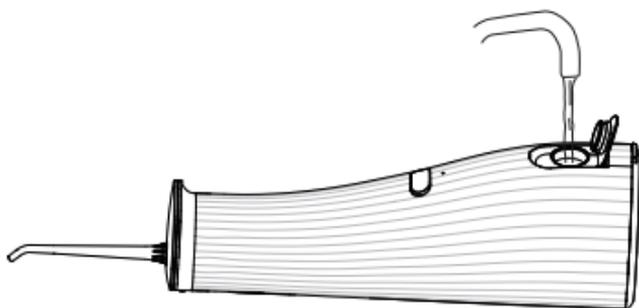05 OPERATION

- Before using the water flosser, brush your teeth thoroughly with a toothbrush.
- Place the nozzle tip into your mouth at a 45° angle and press the mode button  to select one of the modes. Press . Holding the nozzle at 90° angle, direct the jet horizontally into the area between the teeth along the line of the gum. Move the nozzle from the back teeth to the front teeth. Do not press the tip too hard against the teeth or gums.
- The water flosser will turn off automatically when the tank is empty. To turn the water flosser off sooner, press the button .
- Drain any remaining water from the tank.

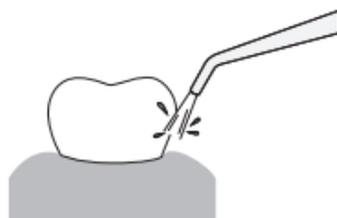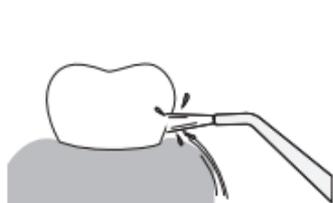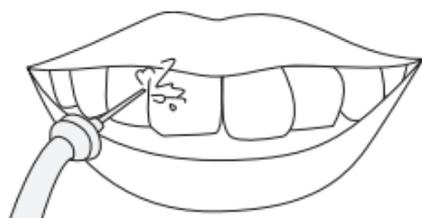
Operating modes

Normal: for gentle daily oral care: recommended for sensitive gums and when the gums have become accustomed to the water pressure in SOFT mode.

Soft: for gentle care with low water pressure: recommended for sensitive gums and in the first few days of irrigator use, until the gums are used to the water pressure.

Pulse: for spot cleaning and gums massage.

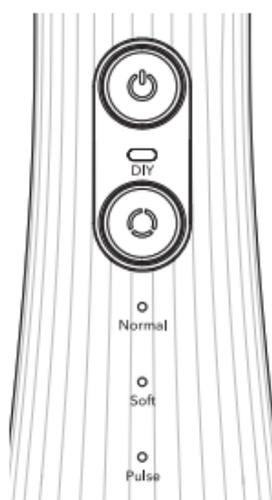DIY (my program): the mode allows you to adjust the intensity of the water supply. A total of 7 intensity levels are available. Press the mode button to choose the DIY mode and hold the button to adjust the intensity level.

Recommendations for use

- Use the water flosser over the sink, as splashes may fly from the water flosser.
- At the beginning of use, it is recommended to use the water flosser once a day or every other day.
- Slight bleeding of the gums may occur during the first few uses. Bleeding should subside after 1-2 weeks of regular use; if gums continue to bleed, discontinue use and contact your dentist.
- Water flosser is not a substitute for brushing teeth with toothpaste and flossing. It is only effective when used together with these procedures.
- Too frequent use of an irrigator can thin tooth enamel and increase tooth sensitivity.

06 CLEANING AND MAINTENANCE

Wipe the body of the appliance with a clean, damp cloth and wipe dry.

- Rinse the nozzles under cool water and soap and water. It should be done after each use. It is recommended to replace the nozzles every three months.
- Drain the reservoir, then turn the appliance on and let it work out the remaining water.

- In case you plan not to use the appliance for a long time, you should remove the nozzle, drain the water and fully charge the appliance.
- Store the appliance in a cool, dry place out of the reach of children.

WARNING!

- **Do not use abrasive cleaning agents.**
- **Do not immerse the appliance in water!**

07 TROUBLESHOOTING

Problem	Possible cause	Solution
Water flosser does not turn on.	The water flosser is completely discharged.	Charge the water flosser according to these the instruction manual.
Low water pressure.	Low battery level.	Fully charge the water flosser.
	The nozzle is clogged.	Replace the nozzle.
Water is dripping from the tank.	The plug of the water inlet is not closed.	Close the cover tightly.
When plugged in, the water flosser does not charge.	The USB cable is damaged.	Contact an authorized service center.
	The appliance is not connected to the power supply properly.	Ensure that the appliance is connected to the power supply properly.
The battery charge runs out too quickly after being fully charged.	The water flosser is used at too low an ambient temperature.	Use and charge the water flosser at temperatures between 15 and 35 °C.

WARNING!

Never repair the appliance yourself. Always contact an authorized service center.

08 DISPOSAL

To protect the environment, do not dispose of appliance with household waste. Bring the appliance to specialized points for further disposal. For more information on the disposal of this product, contact your local authority, your household waste disposal service, or the retailer.

The battery is not intended to be replaced by the user. The battery contains materials that are hazardous to the environment. The battery must be removed from the appliance before disposal. The instrument must be unplugged when the battery is removed. The battery must be disposed of safely in accordance with the laws of the country of sale.

09 TRANSPORT AND STORAGE

- Store the appliance in a clean, closed, dry room at an ambient temperature between 5 and 40 °C and a relative humidity of no more than 70 %.
- Transportation of the device can be carried out by all means of transport, in accordance with the rules for the carriage of goods. When transporting in the original packaging, protect the appliance from damage, preventing the exposure to precipitation and aggressive environments.
- When transporting the appliance without original packaging, use bubble wrap and corrugated cardboard packaging. Wrap each element of the appliance in a film and place it vertically in the container. The free distance between the elements of the product must be laid with air bubble wrap or folded cardboard.

10 TECHNICAL SPECIFICATIONS

Electric sonic tooth brush

- Model: CRP100
- Power: 10 W
- Battery: Li-ion, 2000 mA·h, 3,7 V
- Input: 5 V = 1 A
- Water pressure: 50-150 PSI
- Pulsation frequency: up to 1600 pulses per minute
- Water tank: 310 ml
- Charging time: 4 h
- Protection rating: IPX7
- Product dimensions: 293 x 85 x 65 mm
- Net weight: 260 g

11 SCOPE OF DELIVERY

- Water flosser
- Nozzles - 3 pcs
- USB cable
- Pouch
- Instruction manual

12 WARRANTY

Organization for warranty service and acceptance of claims:
LLC "MVM", Russia, 105066, Moscow, Nizhnyaya
Krasnoselskaya str., 40/12, building 20, floor 5, room II,
room 3.

Phone: +7 (495) 644-28-48

Information about the warranty conditions:

<https://www.eldorado.ru/spv>

or by phone Eldorado - 8 (800) 250 25 25

<https://www.mvideo.ru/spv>

or by phone M.Video - 8 (800) 600 77 75

The seller's warranty is 12 months from the date of sale
(transfer) of the goods to the Customer.

13 ADDITIONAL INFORMATION

Manufacturer: Iconbeauty (Shenzhen) Holding Ltd.

Address: Room 703, Block 1, Idee Monto Industry Park, 4th Industrial Zone, Shutianpu Community, Matian Street, Guangming District, Shenzhen City, Guangdong Province, China

Made in China

Importer/Manufacturer's authorized company: LLC "MVM", Russia, 105066, Moscow, Nizhnyaya Krasnoselskaya str., 40/12, building 20, floor 5, room II, room 3.

E-mail: cert-log@mvideo.ru

Phone: +7 (495) 644-28-48

Installed by the manufacturer in accordance with paragraph 2 of article 5 of the Law of the Russian Federation "On Consumer Protection" the service life of this product is 3 years from the date of sale, provided that the product is used in accordance with the rules and recommendations set forth in this instruction manual and applicable technical standard.

There are no special conditions for the sale.

The date of manufacture is indicated on the appliance where MM is the month of manufacture and YYYY is the year of manufacture.

The product is certified, meets the requirements:

TR TS 020/2011 "Electromagnetic compatibility of technical means"

TR EAEU 037/2016 "On the restriction of the use of hazardous substances in products of electrical engineering and radio electronics".

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