CARRERA

DESIGNED IN GERMANY

CARRERA

DESIGNED IN GERMANY

EN
WATER
FLOSSER

Instruction manual

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№106

Welcome to the **CARRERA** club!

The story of **CARRERA** starts in the 80' in Germany as a lifestyle across two product categories: personal care and home appliances.

CARRERA engineers and designers were inspired by the legendary Carrera Panamericana, the longest and the most challenging race in the 50'. Only the racers who were ready to go beyond the limitations of their own bodies and materials stood a chance of winning.

Today, that same philosophy drives **CARRERA** to develop unique and powerful products.

Each product is made of hi-tech components and quality materials outlining power, quality, innovative features and effective design.

This approach resonated with millions of people, boosting **CARRERA** brand recognition in Europe up to 90 %. **CARRERA** products have become an embodiment of German precision and craftsmanship.

Today, **CARRERA** continues to invest into research and development. Each **CARRERA** product provides users with a sense of luxury, comfort and efficacy.

CARRERA is not just about tools. It is about a lifestyle.

CONTENTS

| | |
|--------------------------------------|-----------|
| 01 SAFETY PRECAUTIONS | 3 |
| 02 INTENDED USE | 4 |
| 03 APPLIANCE OVERVIEW | 5 |
| 04 BEFORE USE | 6 |
| 05 OPERATION | 8 |
| 06 CLEANING & MAINTENANCE | 11 |
| 07 TROUBLESHOOTING | 12 |
| 08 DISPOSAL | 13 |
| 09 TRANSPORT AND STORAGE | 13 |
| 10 TECHNICAL SPECIFICATIONS | 14 |
| 11 SCOPE OF DELIVERY | 14 |
| 12 WARRANTY | 14 |
| 13 ADDITIONAL INFORMATION | 15 |

01 SAFETY PRECAUTIONS

Read this instruction manual carefully before using the appliance. Keep this instruction manual, the packing box and, if possible, other elements of the appliance's packaging. The manufacturer shall not be liable for any damage or malfunction of the appliance or for any damage resulting from incorrect or improper use of the appliance.



- Check the appliance and USB cable regularly for damage. Do not use the product if any damage is detected.
- Do not use the product for other than its intended use.
- Do not use the product outdoors.
- Do not immerse the product or USB cable in water or other liquids.
- Do not expose the instrument or USB cable to high temperatures, direct sunlight, dampness, or other adverse conditions.
- If you have recently undergone oral surgery or are undergoing oral treatment, it is not recommended to use the appliance without prior consultation with a doctor.
- If your gums start bleeding during use, you should stop using the product until you consult your doctor.
- Remove jewelry (rhinestones, grills, etc.) from the teeth before using the appliance.
- Always turn the appliance off before changing the nozzle or filling the water tank.
- Do not turn the appliance on when the water tank is empty. Do not fill the water tank with hot water.
- Do not leave the appliance running unattended, turn it off immediately after use.
- Be sure to pour out any remaining water after use. Clean the appliance after each use.
- Do not direct the water stream into the eyes, under the tongue, into the ear, nose or other sensitive areas.
- Do not use accessories that are not included in the scope of delivery.
- Do not subject the appliance to mechanical impact. Do not apply excessive pressure to the nozzle. Do not use deformed or mechanically damaged nozzles.
- Do not allow children to play with the appliance, the USB cable or the packaging.

- The appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience or knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety.
- Do not attempt to repair the appliance or replace any parts yourself. In case of malfunctions, contact the nearest service center.

02 INTENDED USE

The portable water flosser is designed for oral care and cleaning of the interdental space. The appliance is intended for use in commercial or domestic areas in accordance with the intended use and the safety precautions described in these operating instructions.

03 APPLIANCE OVERVIEW

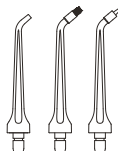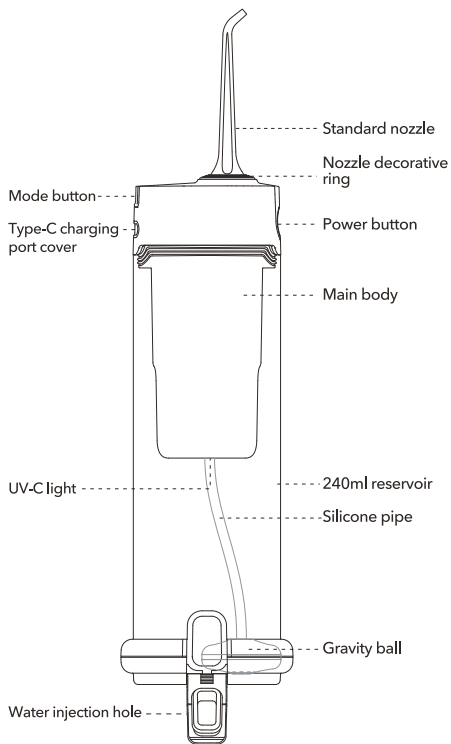3 Nozzles


USB cable

04 BEFORE USE

Charging the appliance

- Charge the water flosser for 3 hours before using it for the first time. This charging is necessary every time the water flosser has not been used for a long time.
- For charging remove the silicone plug on the body and connect the water flosser to the power supply via USB cable.
- During use, charge level is determined by the color of the battery indicator on the display : if the charge level is below 10 %, the indicator will blink red. If the charge level is above 10 %, the indicator light will be steady.
- During charging, the battery indicator will flash yellow until the battery level reaches 100 %. One full battery charge is enough for 100 days of daily use.
- It is not possible to use the water flosser while it is charging.

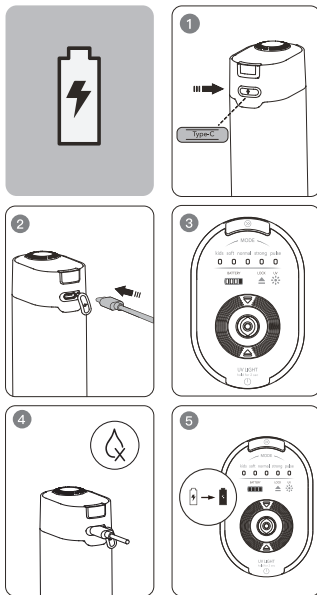
WARNING!

- To maximize battery life, do not leave the battery discharged for long periods of time.
- Do not touch the USB cable with wet hands!
- If the appliance or USB cable become very hot during charging, unplug the appliance immediately and contact an authorized service center.


Filling the water tank

- Before first use, rinse the water tank under running water and drain the remaining water.
- Fill the tank with clean water. This can be done in two ways: through the lid in the appliance body or through the lower opening in the tank.
- Close the cover of the opening until it clicks into place.

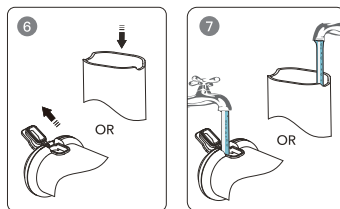Choosing and installing the nozzle

- The set includes 3 interchangeable nozzles:


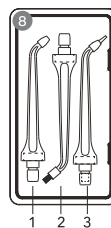
- 1 Standard nozzle for daily oral care.
- 2 Periodontal nozzle for inflamed gums (suitable for periodontitis, periodontal disease and other gum diseases).
- 3 Orthodontic nozzle for cleaning braces as well as bridges, implants and hard-to-reach areas of the mouth.

- Insert the nozzle into the special hole in the top of the housing until it clicks into place.





Make sure it is securely locked in place. If it's not secure in place,  will blink yellow. The water flosser won't start working.

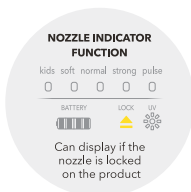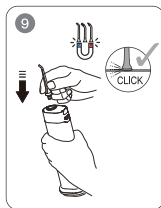- To remove the nozzle, pull it out.






05 OPERATION

- Before using the water flosser, brush your teeth thoroughly with a toothbrush.
- Place the nozzle tip into your mouth at a 45° angle and press the mode power button  to select one of the five modes. Press the power button . Holding the nozzle at 90° angle, direct the jet horizontally into the area between the teeth along the edge of the gum. Move the nozzle from the back teeth to the front teeth. Do not press the nozzle too hard against the teeth or gums.
- The water flosser will run for 2 minutes and then turn off automatically. To turn the water flosser off sooner, press the button .
- Drain any remaining water from the tank.

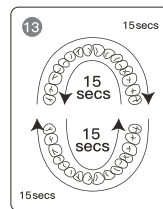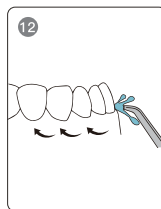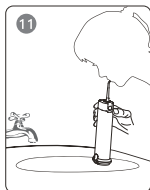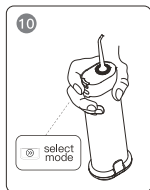Operating modes

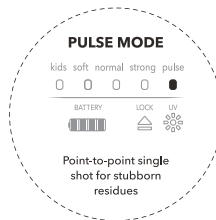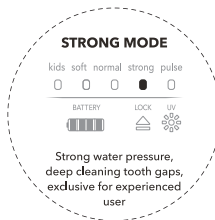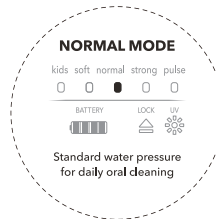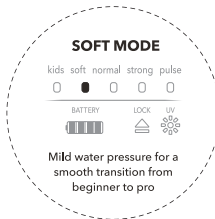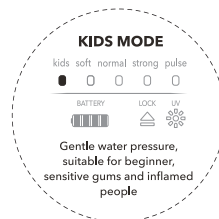
Memory function


The water flosser remembers the last mode used and will be activated automatically the next time you turn the water flosser on. The next time it is switched on, the mode will be activated automatically.

UV-treatment

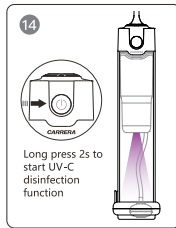UV-treatment allows you to disinfect not only the water tank, but also the nozzles when placed in the water tank. It is recommended to use the function at least once a week.

Hold down the  button for a few seconds to activate the function. The yellow indicator will light up. The function can be activated with the water flosser on or off. Hold down the button  again to deactivate the function. After 5 minutes, the function will be deactivated automatically.


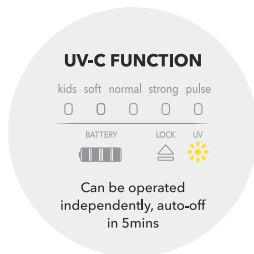
Notes: The UV treatment will be deactivated when separate the tank from the device for safer use.



Recommendations for use

- Use the water flosser over the sink, as splashes may fly from the water flosser.
- At the beginning of use, it is recommended to use the water flosser once a day or every other day.
- Water flosser is not a substitute for brushing teeth with toothpaste and flossing. It is only effective when used together with these procedures.
- Too frequent use of an irrigator can thin tooth enamel and increase tooth sensitivity.
- To prevent scale on the inner surface of the tank, it is recommended to rinse it thoroughly after each use.

06 CLEANING & MAINTENANCE

- Wipe the body of the appliance with a clean, damp cloth.
- Rinse the nozzles under cool water and soap and water - this should be done after each use. The nozzles can be stored in the water tank. It is recommended to replace the nozzles every three months.
- If necessary, clean sediment or scale in the water tank with a brush.
- During storage, it is recommended to keep the lid of the water tank open to prevent unpleasant odor.
- If you do not plan to use the appliance for a long period of time, remove the nozzle. Fully charge the appliance.
- Store the appliance in a cool, dry place out of the reach of children.

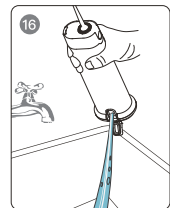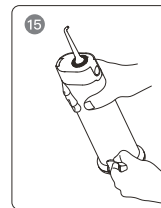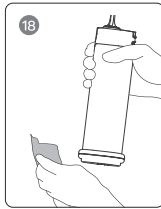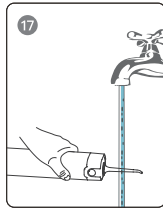WARNING!

**Do not use abrasive cleaning agents.
Do not immerse the appliance in water!**

07 TROUBLESHOOTING

| Problem | Possible cause | Solution |
|--|--|---|
| Water flosser does not turn on. | The water flosser is completely discharged. | Charge the water flosser according to these the instruction manual. |
| Low water pressure. | Low battery level. | Fully charge the water flosser. |
| | The nozzle is dogged. | Replace the nozzle. |
| | The gravity ball filter is clogged. | Clean the filter. |
| Water is dripping from the tank. | The plug of the water inlet is not closed. | Close the cover tightly. |
| | The sealing ring of the water inlet cover is damaged. | Contact an authorized service center. |
| When plugged in, the water flosser does not charge. | The USB cable is damaged. | Contact an authorized service center. |
| | The USB cable is connected incorrectly. | Check the connection. |
| The battery charge runs out too quickly after being fully charged. | The water flosser is used at too low an ambient temperature. | Use and charge the water flosser at temperatures between 15 and 35 °C. |
| The irrigator has switched off during operation and the yellow light on the panel is flashing. | The overvoltage protection has been activated. | Wait 15 seconds. When the light goes out, turn the water flosser back on and continue using it. |

WARNING!

Never repair the appliance yourself. Always contact an authorized service center.

08 DISPOSAL

To protect the environment, do not dispose of appliance and the batteries with household waste. Bring the appliance and batteries to specialized points for further disposal. For more information on the disposal of this product, contact your local authority, your household waste disposal service, or the retailer.

The battery is not intended to be replaced by the user. The battery contains materials that are hazardous to the environment. The battery must be removed from the appliance before disposal. The appliance must be unplugged when the battery is removed. The battery must be disposed of safely, in accordance with the legislation of the country of sale.

09 TRANSPORT AND STORAGE

- Store the appliance in a clean, closed, dry room at an ambient temperature between 5 and 40 °C and a relative humidity of no more than 70 %.
- Transportation of the device can be carried out by all means of transport, in accordance with the rules for the carriage of goods. When transporting in the original packaging, protect the appliance from damage, preventing the exposure to precipitation and aggressive environments.
- When transporting the appliance without original packaging, use bubble wrap and corrugated cardboard packaging. Wrap each element of the appliance in a film and place it vertically in the container. The free distance between the elements of the product must be laid with air bubble wrap or folded cardboard.

10 TECHNICAL SPECIFICATIONS

- Model: CRP106
- Battery: Li-ion, 2500 mA·h, 3,7 V=
- Input: 5 V = 2 A
- Water pressure: 276-896 kPa
- Pulsation frequency: 1100-1600 pulses per minute
- Water tank: 240 ml
- Charging time: 3 h
- Protection rating: IPX7
- Product dimensions: 305 x 77 x 54 mm
- Net weight: 305 g

11 SCOPE OF DELIVERY

- Water flosser
- Nozzles - 3 pcs
- USB cable
- Pouch
- Instruction manual

12 WARRANTY

Organization for warranty service and acceptance of claims: LLC "MVM", Russia, 105066, Moscow, Nizhnaya Krasnoselskaya str., 40/12, building 20, floor 5, room II, room 3.
Phone: +7 (495) 644-28-48

Information about the warranty conditions:
<https://www.eldorado.ru/spv>
or by phone Eldorado - 8 (800) 250 25 25
<https://www.mvideo.ru/spv>
or by phone M.Video - 8 (800) 600 77 75

The seller's warranty is 12 months from the date of sale (transfer) of the goods to the Customer.

13 ADDITIONAL INFORMATION

Manufacturer: Guangdong Energetic Home Appliance Co., Ltd

Address: No.20 Changjiang Road, Sanshan New Town, Guicheng, Nanhai District, Foshan City, Guangdong, China

Made in China

Importer/Manufacturer's authorized company:
LLC "MVM", Russia, 105066, Moscow, Nizhnaya Krasnoselskaya str., 40/12, building 20, floor 5, room II, room 3.
Phone: +7 (495) 644-28-48
E-mail: cert-log@mvideo.ru

Installed by the manufacturer in accordance with paragraph 2 of article 5 of the Law of the Russian Federation "On Consumer Protection" the service life of this product is 3 years from the date of sale, provided that the product is used in accordance with the rules and recommendations set forth in this instruction manual and applicable technical standard.

There are no special conditions for the sale.

The date of manufacture is indicated on the appliance where MM is the month of manufacture and YYYY is the year of manufacture.

The product is certified, meets the requirements: TR TS 020/2011 «Electromagnetic compatibility of technical means»
TR EAEU 037/2016 «On the restriction of the use of hazardous substances in products of electrical engineering and radio electronics».

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